



## 70PL Fractional IPL

#### Tentech Item





**10 THERMA**—
Monopolar RF



10 THERA
\_\_\_\_\_
2 Line HIFU



10PL — 168 Cell Fractional IPL



3cm x 3cm CO2 Fractional Laser



Dual Mode
Ultrasound System



TENHI
—
Fractional & Thermal RF



10 50NO





LIPO THERA

Dual Mode

Ultrasound System



High Intensity &
Multiple Light Therapy
Device



PhotoBioModulation
Device



Regevan® LED Mask Dr. Oracle LED Mask 4,474 LED Lights



Mono Thread
[Registered for NMPA]



10 cog



Rhino fill
Premium
Regevan fill
Volume PRO
Regevan fill
Volume UP
Regevan fill
Fine
Monophasic Filler
[for Export]



Regevan Polynucleotide Sodium



Regevan°C



Regevan<sup>®</sup>5
Growth Factor Spray



Regevan®
\_\_\_\_ ALO CARE
Growth Factor Ampoule



### 10PL

#### Fractional IPL

- Blemish · Freckle · Age spot
- Pigmentation · Cafe au Lait · Acne
- Partial lentiginosis Flushing










10PL is an IPL device that significantly lowers the risk of burns by applying the fractional treatment concept to the existing IPL. Also, the tip size is 45mm X 18mm, which allows for quick treatment of large areas





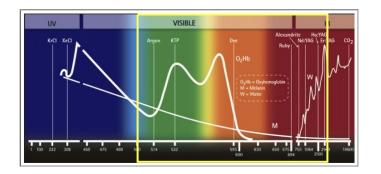

IPL is...


Intense Pulsed Light emits light of multi-wavelength in the form of strong pulse, From 400~1,200nm wavelength and is effective for treating Acne, Pigmentation, Flushing, etc.

#### IPL and Laser Comparison

	IPL	LASER
Wavelength	Multiple wavelength	Single wavelength
Treatment areas/effects	Simultaneously treats lesions such as pigmentation, vascular, pore, elasticity	Treats single lesion
Down time	Daily life activities immediately after treatment	1-10 days depending on the laser
Pain	Anesthetic cream	Anesthetic cream or local anesthesia





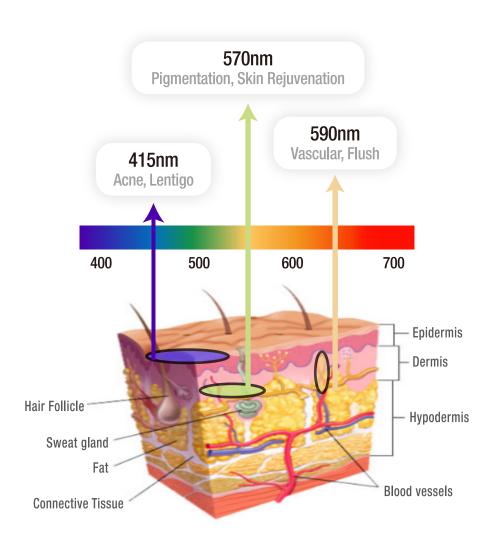





Technical Specification		
Light Source	Intense Pulsed Light (Xenon Lamp)	
Standard Spectra (Wavelength)	415nm / 570nm / 590nm	
Energy Density (Fluence)	$5-35 J/cm^2$	
Spot Size	45 x 18 mm	
Pulse Sequence	Single , Double , Triple pulses	
Pulse Duration	2.0 ~ 35ms	
Pulse Delay	3.0 ~ 60ms	
Repetition Rate	0.5 ~ 1 Hz ± 10% (Max, Energy Fluence)	
Light Conduction	Sapphire Contact Coupling	


System Specification		
Electrical Power	220 V, Single Phase, 50/60 Hz, 16 Amps at 220V AC	
Dimension	378 mm(W) x 560 mm(L) x 1070 mm(H)	
Weight	73 Kg	





#### 1) Fractional IPL allows increased safety

Through 168 Cells (each1.4 x 1.59mm), thermal damage transmitted to the epidermis is minimized, and safe and effective treatment is possible even with higher energy.

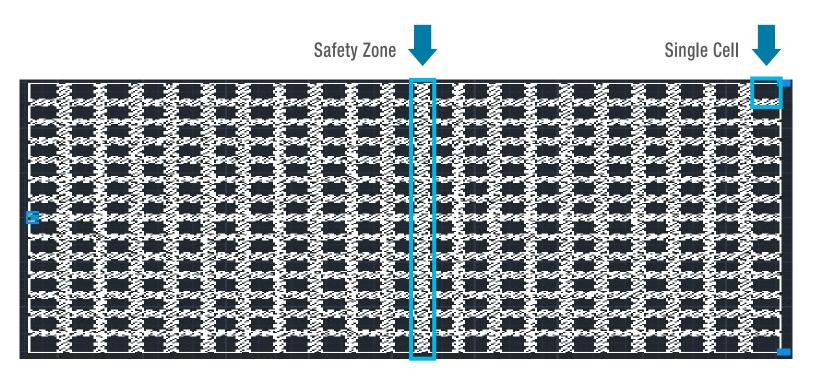





# Conventional IPL 10PL MANAGEMENT





#### 2) Select Wavelength using Turn Dial

Treat various lesions with three wavelengths without exchanging filters!





Wavelength (nm)	Rotary switch symbol
415	4
570	5
590	9





3) 45mm x 18mm tip size reduces treatment time

38x10mm

35x15mm

45x18mm

**Other** 

10PL

4)Fast Procedure with Beam Irradiation Speed of 1Hz or less

5) Treat Blemish using 35J High Power Fractional Beam



#### 10PL is for patients who...





concerned about freckles, age spots,
 blemishes, black spots, and pigmentation

 worried about cafe au lait, nevus spilus, and Partial lentiginosis troubled with redness and acne



### 10PL

All information in this document is subject to copyrights and other intellectual property rights of Tentech Inc. and other licensees. No modification, reproduction or copy of any or any part of this material may be made without prior written permission from Tentech Inc. and licensees.

Additional information is available upon request. The following are trademarks or registered trademarks of Tentech, Inc. in Korea and other countries.

